




Goldfern
MOBILITY™

ZERO GRAVITY LIFT CHAIR

INSTRUCTION MANUAL

Read this manual carefully before using this product

ZERO GRAVITY LIFT CHAIR

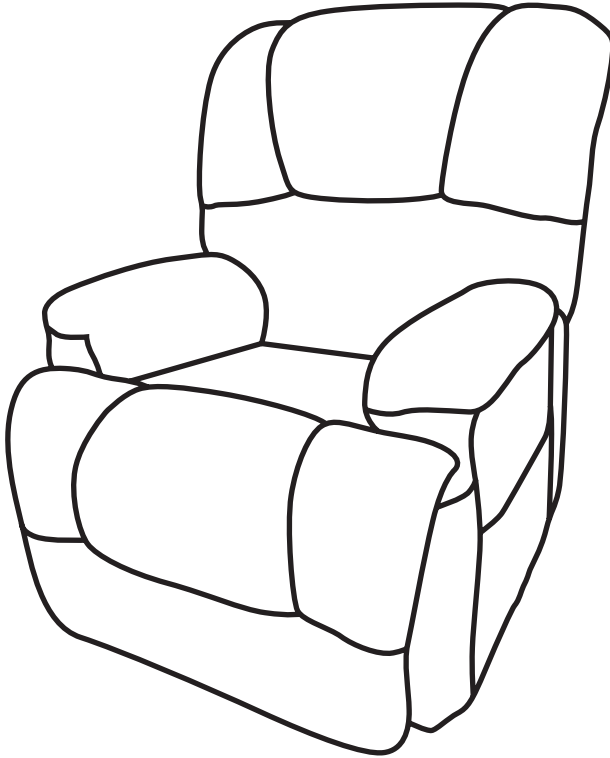

Specifications

Weight

- | | |
|--------------------------|--------|
| • Weight of chair | 57kgs |
| • Maximum loading weight | 150kgs |

Power

- | | |
|-----------------|----------|
| • Rated Voltage | 230-240V |
| • Frequency | 50Hz |
| • Rated Power | 50W |

Chair

- | | |
|--------------------------|-----|
| • Low voltage for safety | 24V |
|--------------------------|-----|

Features

- Comfortable design
- Easy to use
- A complete and easy to use remote control
- Dual motors
- Lift and recline function
- External battery backup

Goldfern Zero Gravity Lift Chair

MOBILITY™

This is an electrical lift chair designed to help people rise with ease from the chair. Before you start using the chair, please read our instruction manual carefully and follow the safety procedures recommended. Keep this manual safe for future reference.

Safety Guidelines

- The chair must be assembled correctly according to the Assembly instruction enclosed.
- The manufacturers of the chair, their agents and distributors are not liable for any damages caused by improper assembly of the chair.
- The chair must be assembled by using all the parts provided.
- The chair must be placed on a flat floor level.
- The chair must be in a dry room.
- Do not expose this chair to direct sunlight, stoves, heaters or any other heating appliances.
- Do not use the chair outdoors.
- When the chair is in storage, disconnect the power plug from the main supply to avoid overheating.
- Do not turn off chair over night.
- Do not use this chair in combination with other therapeutic or beauty equipment such as massage mat, slimming machine or electric blanket.
- Do not allow children or pets to play on or around the chair. This chair is not for children.
- Do not place any heavy object including children on the armrest, footrest, or backrest.
- Do not jump or play on the chair.
- Do not rock the chair, it is not a rocking chair.
- Do not use the chair when your body is wet or damp.
- Do not sleep on the chair as opposed to a bed.
- Do not connect or disconnect the power plug to/ from the electrical outlet.
- When disconnecting the plug, pull it out by holding the plug and not by the cord.
- Do not connect or disconnect the plug from the electrical outlet with wet hands. It may result in electrical shock or machine trouble.
- To avoid electrical shock, the power plug and cord must always be dry. Never wet the power plug and cord or touch them with wet hands.
- Do not place the power cord under any objects including the chair as it may damage the cord and result in fire or electrical shock.
- Do not let the mains power cord hang down from the chair. Must be kept clear of all moving parts.
- Do not operate the remote control with wet hands or with sharp or hard objects.
- Before using the chair, check the power cord or plug is not damaged. Do not operate the chair if the power cord or plug is damaged. Do not try to change the damaged plug or cord by yourself. Always contact a qualified technician or the After Sales Service for inspection and repair.
- Do not sit, step or place heavy objects on the remote control.
- Do not operate the remote control with excessive force.
- Do not operate the chair with damaged or broken remote control. It may result in electrical shock. The remote control must be replaced by a new one and must be done by a qualified technician.
- Do not sit or lie on your side or stomach during operation of the chair. You must always sit or lie upright and on your back.
- Do not put your fingers or any object between the seat and footrest.
- Before cleaning or repairing the chair, make sure the main switch is turned off and the power plug disconnected from the electrical outlet.
- Do not use this chair in humid or dusty environment, as it may cause electrical shock or machine trouble. The chair must be in a dry room.
- Do not use this chair in room temperature higher than 40°C.
- Do not wash the chair with water.
- Do not use the chair if the upholstery or any part of the chair is damaged.
- The chair can bear a maximum weight of 150kg. Person(s) above this weight should not sit on this chair as it may damage the chair.
- Do not remove the upholstery or the back cover from the chair. Touching the inside mechanisms of the chair may result in electric shock.
- The manufacturer of this chair, their agents and distributors are not liable for any damages caused to person(s), object(s), animal(s) or property(s) resulting from improper and misuse of the chair.
- Keep this instruction book for future reference.

When to stop using the chair:

- If water is spilt on the chair. Continuous usage may result in electric shocks.
- If the upholstery or any part of the chair is torn or damaged.
- If the chair has a malfunction during operation.
- If there is a power failure. Injury may occur if the power is restored unexpectedly.
- If you suffer from abnormal or severe pain during the operation of the chair.
- In any of the above events switch OFF immediately and disconnect the power cord from the power outlet.


Further notes:

- Check for tears in the upholstery cloth regularly.
- Check all wires, plug and the remote control regularly for any damage.
- Check all wires are free from obstruction and will not be crushed when lifting or lowering the chair.
- Power transformers are all independently tested and tagged in New Zealand, to New Zealand standards.

Assembly Instructions

Unpackaging the Lift Chair

Cut open the top of the box, then carefully flip the box upside down and lift box up to remove.

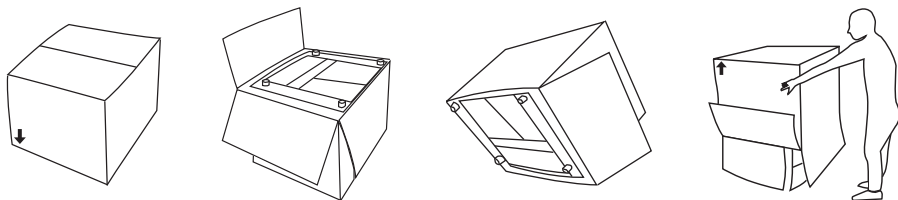


! IMPORTANT

Do not attempt to lift the Lift Chair out of the box. Lifting the chair may cause injury due to its weight, or may damage the chair itself.

Assembling the Chair

The Lift Chair is packed partially assembled, with separate backrest and chair base parts.


Step 1: Slide the backrest down to engage the sleeve with the base bracket, ensuring it locks into place.

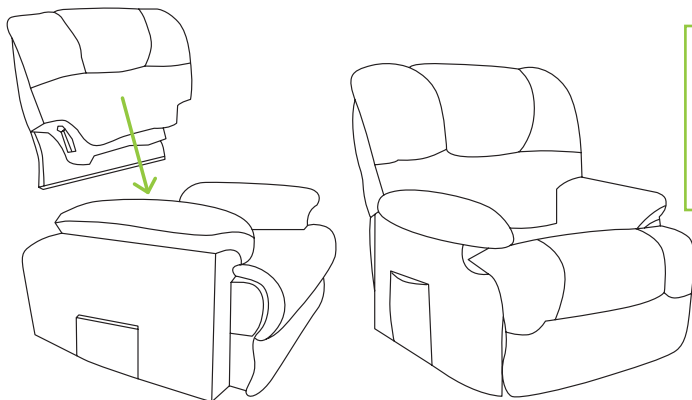

! IMPORTANT

Never attempt to pick up or move the chair by holding on to the chair cushions or arms. This may stretch or tear the material of the chair.

Step 2: Connect cables from the backrest to the wire from the seat as shown below.

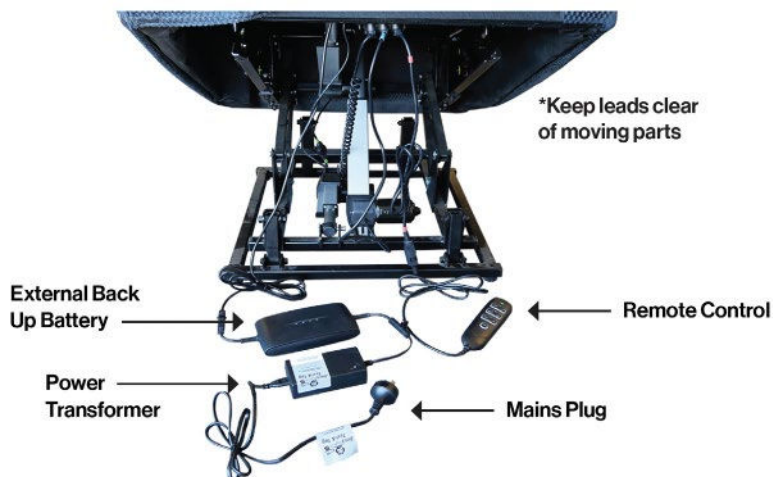

! IMPORTANT

Ensure that there are no obstructions between the chair base and footrest when using the lift function. Before using the chair, ensure all cords are free of any moving parts.

ALWAYS place the chair on a flat and smooth floor. **DO NOT** place the chair on an uneven surface.

User Instructions

Positioning the Chair

Ensure that the chair is positioned as described to allow safe reclining and lifting movements.



! CAUTION

- Make sure there is no obstruction behind, in front/ or above the chair for a safe recline and lift during operation.
- No finger or object should be placed between the chair and the footrest.
- Before using the chair, ensure cords are free of any moving parts.
- Once the chair is assembled, place in a dry room, with sufficient space in front and behind it.
- Make sure there is an allowance of more than 70cm behind the chair and more than 40cm in front.
- Place the chair on a flat and smooth floor. Do not place the chair on uneven floor.

Remote Control

The remote control can be set up for left or right-handed use.



To sit in the chair

1. While standing in front of the chair, press and hold the **Mid Right** button to raise the chair to a suitable height for you to sit down.
2. Press and hold the **HOME** button to fully lower the chair. Stop pressing the button once you reach a comfortable position - the chair will continue to transition from assisted standing to a neutral sitting position.

Footrest:

Raise - press and hold the **Mid Left** button / Lower - press and hold the **Mid Right** button.

Backrest:

Recline - press and hold the **Top Left** button / Raise Upright - press and hold the **Top Right** button.

Zero Gravity:

Tilt Backwards - press and hold the **Bottom Left** button / Tilt Upright - press and hold the **Bottom Right** button.

To get out of the chair

1. Return chair to a neutral seated position by holding **HOME** button.
2. Press and hold **Mid Right** button to raise the chair to a suitable height for you to stand easily.

Zero Gravity



! IMPORTANT

For safety purposes, chair should be left in the neutral sitting position (footrest lowered, backrest upright) when not in use. To prevent children, pets or obstructions from becoming trapped underneath the chair.

Repairs & Maintenance

As a family-owned and operated New Zealand company, we are committed to delivering quality products that improve daily living with reliability and care. For your safety and the best performance and longevity of this product, please follow the instructions below when your lift chair requires repairs or maintenance.

The symbols below are used throughout this owner's manual to identify warnings and important information. It is very important for you to read and understand them completely.

WARNING

Indicates a hazardous situation that could lead to serious injury or death if it is not avoided.

! IMPORTANT

Indicates a hazardous situation that could result in damage to property if it is not avoided.

CAUTION

Indicates a hazardous situation that could result in minor or slight injury if it is not avoided.

Product Repairs

If your chair requires repairs, please contact the After Sales technician to arrange for a professional to carry out the repairs/servicing required. Before any repairs are carried out, always ensure that the power is turned off and the power plug is removed from the electrical outlet.

- The chair must be serviced/repaired by a professional.
- Do not attempt to repair any parts of the chair by yourself.
- Do not change the power plug by yourself.
- Do not disassemble the chair or remove any inside mechanisms by yourself.

Product Care & Maintenance

Before cleaning your chair, always ensure that the power is turned off and the power plug is removed from the electrical outlet.

- The upholstery of the chair must be cleaned with a dry cloth only.
- Do not use any cleaning agents, thinners, or bleaching agents to clean the upholstery.
- Do not wash the upholstery.
- Ironing of the upholstery is strictly prohibited.
- The handles and chair must be cleaned with a dry cloth. Never use water to wash/clean any part of the chair.

CAUTION

Do not remove the upholstery for cleaning. Touching inside chair mechanism of the chair may result in electrical shock.

! IMPORTANT

The remote control **must only** be cleaned with a dry cloth. Never use water or a wet cloth, as this may cause damage to the electrical components inside the remote control.

Product Storage

Before storing your chair, ensure the power is turned off and the power plug is removed from the electrical outlet.

- Ensure that the chair is stored in a clean state. Wipe off any dust and dirt from the chair surface before storing.
- Do not store the chair in high-humidity environments.
- If the chair is not in use for an extended period of time, cover it with a clean sheet to protect it from dust.
- Do not store the chair in direct sunlight or high-temperature environments as the fabric colour may fade.

NOTES:



Our team understands quality rehabilitation equipment and we can provide you with pieces that will become an integral and trusted part of your daily routine. We are a family business dedicated to providing you with equipment you deserve at a price that makes sense.

The Goldfern Mobility story began as we walked along our production premises and learned we can supply quality wheelchairs for a reasonable cost. The more we researched equipment needs in the disability sector, studied the rising costs in the market, and listened to client stories, the more passionate we became.

Our goal is to be affordable - but without compromising the functionality or safety of the product. Safety is paramount. With this in mind, we have taken it upon ourselves to lead the way in NZ and, where possible, ensure our products are tested and approved to internationally recognised safety standards. Until such time as these standards become compulsory by law, Goldfern Mobility will continue to invest in testing standards.

Over the years we have listened and learned from our clients. We are intentionally selective with our stock and our size to ensure we meet an expected standard of excellence. We work to modify our equipment as necessary so everything meets our 3-point standard: Safe. Functional. Affordable. We want to do the same for you.

No matter where you are in the country, you can locate your nearest Goldfern Mobility retailer online. We are eager to hear from you and help meet your needs by providing superior quality products for a reasonable price.

We look forward to hearing from you.

hello@goldfernmobility.co.nz | www.goldfernmobility.co.nz | 0800 100 971



208 John F Kennedy Drive, Milson,
Palmerston North 4414
New Zealand

www.goldfernmobility.co.nz

hello@goldfernmobility.co.nz | 0800 100 971

