

ECO-LITE ROLLATOR

INSTRUCTION MANUAL



Read this manual carefully before using this product



Safety Notice & Caution

- Consult your doctor, therapist or other qualified professional to determine if this rollator will suit your particular needs.
- For proper positioning and use of your Rollator, imagine a line drawn across the back of the Rear Wheels. Part of one foot should always be slightly ahead of that line as you walk.
- This rollator is for individual use only.
- Rollator is designed for indoor & outdoor use.
- Regularly check the brakes for proper operation.
Pulling up on the brake handle will slow the rollator.
Pushing down on the brake handle will lock the wheels.
Pull up on the brake handle to release the brake lock.
- Do NOT use as a wheelchair or transport chair.
- Rollator must be fully opened and in a locked position when being used as a walking aid.
- Wear proper shoes while using this rollator.
- Do NOT sit on the rollator and push yourself around using your feet to propel.
- Do NOT stand on seat.
- Before sitting on seat, be sure the rollator is fully opened, and the brake is in the locked position.
- Do NOT use on stairs or escalators.
- DO NOT sit on the rollator if the ground is sloped or uneven.
- Exercise caution when walking on sloped surfaces.
- Do NOT use if the rollator is not functioning properly.
- Maximum user weight 126kg
- When storing/transporting your rollator do not place anything on the rollator, as damage may occur.
- Make sure your fingers are clear of the folding or pinch point, whenever the product is folded or unfolded.

Set Up Instructions

! IMPORTANT

Ensure that all bolts are securely fastened before using this rollator! Brakes must be in the parked/locked position before sitting on the rollator.

Set the handlebars at the required height

To determine this, stand behind the rollator with shoulders relaxed, arm at your side. The hand grips should be set to a height which is level with the crease of your wrist. Consult your physiotherapist or other qualified healthcare professional for assistance to determine the correct height to set the handlebars.

The holes in the frame are hexagonal in shape to accommodate the head of the screw. Place the screw through the aligned holes. Do up the knob ensuring the bolt head is on the inside of the frame.

Check the wheels to ensure they are in working order and roll freely.

The brake cables should be flush along the legs of the rollator. Push the cable into the clip at seat level if not already clipped. Ensure excess brake cable loops out the front of the rollator.

Folding The Rollator


- To fold the rollator, remove the basket (if present). Remove contents of bag (if any).
- If a tray is fitted fold it down the front of the rollator.
- Lift the seat up.
- Release the plastic safety latch on the sidebar of the main frame.
- Pull up on the centre bar from inside the storage bag.

Storage & Transportation

Fold your rollator for easier storage and transportation.

Do not rest anything on your folded rollator when being stored or transported as damage may occur.

Usage

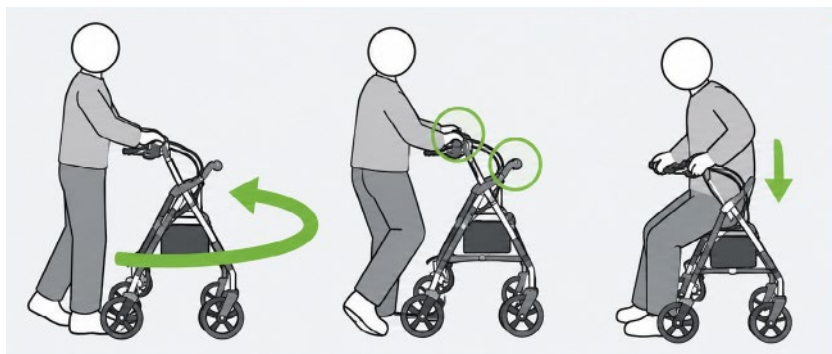


- Do not support your weight on the handle attempting to sit.
- Use the handles to stabilize your body as you carefully turn around.
- Keeps your hand on the handles and slowly lower your weight onto the seat.




Do NOT use as a wheelchair or transport chair.


Brake Usage

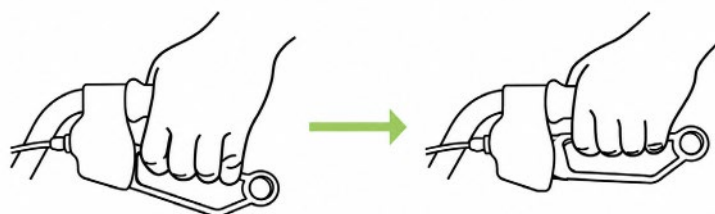

Pull up on brake handle to slow down or to release the parking brake.

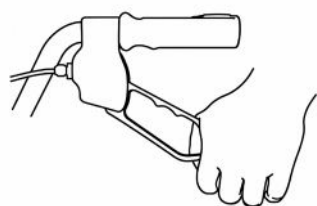

Push down on brake handle to lock the brake into PARK.

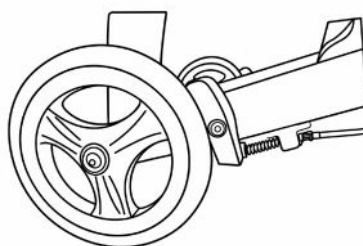

Adjust brake tension, Loosen the locking nut, and turn the adjustment screw clockwise/counter clockwise.

Maintenance

Use a non-abrasive detergent or cleaner with warm water.

Periodically check knobs and screws to ensure they are tight.

Periodically test the brakes by both walking and applying the brakes and by stopping. Push the brakes down into a locked position and then try pushing the walker. If the brakes do not lock the wheels, the brakes need to be adjusted.

Warranty

We warrant this product against defects in material and/or workmanship for a period of one year from date of purchase. Those items not covered under warranty are wear items, such as hand grips, brake cables, and wheels.

We will not warranty the rollator if it has been damaged due to improper installation, modification, and/or misuse.

Service/Labour charges are the discretion on the reseller and are not the responsibility of Goldfern Mobility. Our sole liability extends to repair or replacement of parts or components as previously specified.

ISO Approved

This product has been tested and approved to an internationally recognised safety standard by an independent testing laboratory.





208 John F Kennedy Drive, Milson,
Palmerston North 4414
New Zealand

www.goldfernmobility.co.nz

hello@goldfernmobility.co.nz | 0800 100 971

